



MIRION
TECHNOLOGIES



DOSE CALIBRATORS

CRC[®]-55tPET

Dose Calibrator

*For Speed and Accuracy in
Measuring PET Isotopes*

FEATURES

- 8 in. color VGA touch screen display
- USB/PC Communications
- Software upgrade via USB or Flash Drive
- USB printer capability
- Single or optional second plug and play chamber with mixed chamber capability
- Chamber can be placed 100 feet from the readout unit
- Bilingual (English/Spanish) Software
- On-screen display of Nuclide, Activity, Unit of Measure and Calibration Number
- Large character, high visibility display
- Full alpha numeric touchpad
- Built-in dose calibration, quality control and self diagnostics
- Has a maximum activity up to 20 Ci of F-18
- Compatible with Nuclear Medicine Management Systems via USB
- Optional printer for regulatory records and patient labels for syringes and vials
- Over 80 nuclide symbols and half-lives in memory
- Automated Geometry and Linearity Testing

DESCRIPTION

The new Capintec™ CRC-55tPET dose calibrator combines the ultimate features with the quality and versatility that you've come to expect in one of the industry's finest packages. Using the menu-driven, color touch screen interface, the CRC-55tPET calibrator's reduced chamber pressure and increased bias voltage increases the maximum activity range for high energy PET isotopes. The CRC-55tPET unit's performance and reliability has been designed to meet the demanding needs for speed and accuracy in the preparation and measurement of doses in any laboratory environment.

The enhanced 55t software supports dual chamber technology with an innovative functional design that includes a large, easy-to-read display that specifies Nuclide Name, Number, Activity, Unit of Measure and Calibration Number.

Entering data through the custom touch screen interface is fast and includes 28 programmable keys. The user can select from 80 nuclides by simply selecting the nuclide symbol on the touch screen interface.

Additional capabilities include storage of reference sources in memory that automatically decay corrects for today's time and date. Automated quality control tests and self-diagnostics are built-in with automatic zero and background subtraction making the CRC-55tPET calibrator exceptionally easy to use. An optional printer enables the CRC-55tPET unit to print full size records and patient tickets with peel off labels for vial and syringe identification.

Se habla español? New to the 55t family is bilingual English/Spanish software!

Features including USB/PC communications, printer capabilities, USB flash drive software for upgrades and plug-and-play chamber make the CRC-55tPET calibrator an integral tool in improving your department's efficiency.

Count on us for excellence in energy measurement and unsurpassed customer service, training and support.

CRC-55tPET | DOSE CALIBRATOR

OPTIONAL COMPONENTS




CRC-PS Positron Shield



Epson Roll Printer




Epson Ticket Printer

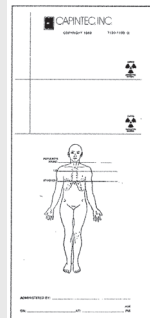

CRP-200 Dose Tickets & Labels



Environmental Shield

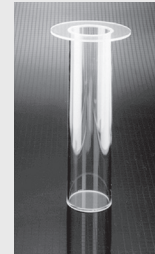

Liner



Dipper

SPECIFICATIONS

Physical:

Console Dimensions

- Height: 24 cm (9.5 in.)
- Width: 23 cm (9.0 in.)
- Length: 27 cm (10.5 in.)
- Weight: 3.4 kg (7.5 lb)

Chamber Dimensions

- Height: 43.8 cm (17.25 in.)
- Diameter: 17.2 cm (6.76 in.)
- Weight: 13.6 kg (30 lb)
- Well Diameter: 6.1 cm (2.4 in.)
- Well Depth: 25.4 cm (10.0 in.)
- Cable Length¹: 3.7 m (12 ft)

Cables

- Power¹: 1.8 m (6 ft)
- Printer: 1.8 m (6 ft)

Ionization Chamber

- Type: Thin wall, deep well
- Fill Gas: 5 atm Ultra-Pure Argon

Measurement Range

- Type: Auto Ranging
- Activity: 20 Curies of F-18
- Resolution: 0.1 μ Ci (0.01 MBq), max.

Display Screen

- Type: 8 in. VGA LCD Color Touch Screen Display
- Format: Direct reading in Ci or Bq
- Bq/Ci Reading: User selectable or fixed
- Values Displayed: Nuclide name (Atomic symbol, Mass number), calibration number

Electrometer

- Accuracy: Better than $\pm 2\%$
- Linearity: Within $\pm 2\%$
- Response Time: Within 2 sec, 4 to 16 sec for very low activity samples
- Bias Voltage: +500V

¹: Longer cables are available. Consult factory.

Repeatability of Measurement

- Within $\pm 1\%$ within 24 hours during which time the calibrator is continuously in operation

Overall Accuracy

- Accuracy Determined By:
 1. Calibration for the specific nuclide and the sample configuration,
 2. Accuracies of standard sources used for calibration of electrometer

Tests

- Diagnostics: Full test of program, system memories
- Daily: Auto Zero, Auto Background, Voltage Test, Data Check, Accuracy and Constancy
- Enhanced: Linearity, Geometry

Nuclear Data

- Nuclide Keys: 28 programmable keys
- System Memory: Over 80 nuclides (Cal number and half-life)

Standard Source Data

- System Memory: Co-57, Co-60, Ba-133, Cs-137, Na-22 standard sources

PC Port

- Interface: RS-232 and USB
- Compatibility: Standard Nuclear Medicine Management Systems

Printer (Optional)

- Interface: RS-232 and USB
- Type: Epson Roll, Epson Slip or Okidata full size dot matrix
- Printing Options: Full size test reports, measured results on tickets

Power Requirements

- 100-240 V ac (50/60 Hz) 100 mA

ORDERING INFORMATION

- CRC-55tPET Dose Calibrator: Item #5130-3235
- CRC-55tPET Multiple Chamber Option: Item #5130-2238
- CRC-PS Positron Shield: Item #7300-2903
- Epson Roll Printer: Item #5430-0058
- Epson Ticket Printer: Item #5430-0100
- CRP-200 Dose Tickets & Labels: Item #7120-1199
- Environmental Shield: Item #7300-2450
- Liner: Item #7300-2004
- Dipper: Item #7300-2005



Mirion Technologies (Capintec), Inc. • 7 Vreeland Road, Florham Park, NJ 07932 USA

OPS-1431 – 01/2020

Copyright © 2020 Mirion Technologies, Inc. or its affiliates. All rights reserved. Mirion, the Mirion logo, and other trade names of Mirion products listed herein are registered trademarks or trademarks of Mirion Technologies, Inc. or its affiliates in the United States and other countries. Third party trademarks mentioned are the property of their respective owners.

MI CAPINTEC, INC.
Part of Mirion Technologies

www.capintec.com
www.mirion.com